

Cellular Senescence Detection Kit - SPiDER-βGal

Technical Manual

Technical Manual (Japanese version) is available at <http://www.dojindo.co.jp/manual/sg03.pdf>

General Information

DNA damages in normal cells are caused by repeated cell division and oxidative stress. Cellular Senescence, a state of irreversible growth arrest, can be triggered in order to prevent DNA-damaged cells from growing. Senescence-associated β-galactosidase (SA-β-gal), which is overexpressed in senescent cells, has been widely used as a marker of cellular senescence. Although X-gal staining is widely used to detect SA-β-gal, this method has the following disadvantages: 1) requirement of fixed cells due to the poor cell-permeability, 2) low quantitative capability because of the difficulty of discrimination between stained cells and unstained cells, 3) requirement of a long time for staining.

Cellular Senescence Detection Kit - SPiDER-βGal allows to detect SA-β-gal with high sensitivity and ease of use. SPiDER-βGal is a new reagent to detect β-galactosidase which possesses high cell-permeability and high retentivity inside cells. SA-β-gal can be detected specifically not only in living cells but also in fixed cells by using a reagent (Bafilomycin A1) to inhibit endogenous β-galactosidase activity. Therefore, SPiDER-βGal can be applied to quantitative analysis by flow cytometry.

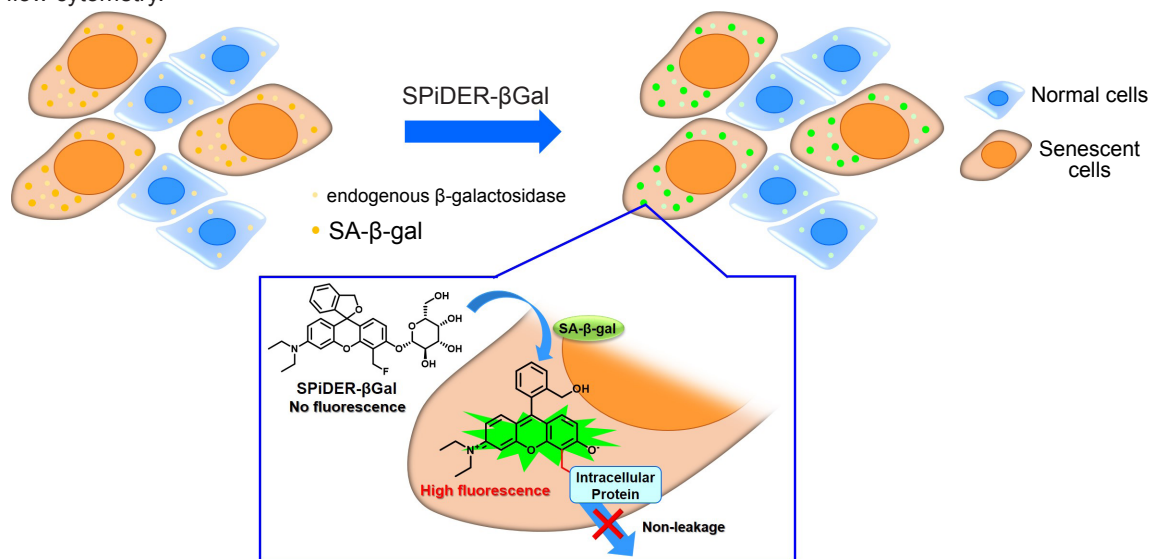


Fig. 1 Detection mechanism of senescent cells by SPiDER-βGal

Kit Contents

10 assays (35 mm dish)

- SPiDER-βGal x 1
- Bafilomycin A1 x 1

※ The Bafilomycin A1 in the tube may be barely visible due to the small amount. Please handle it carefully.

Storage Condition

Store at 0-5°C

Required Equipment and Materials

- Dimethyl sulfoxide (DMSO)
- Culture medium or HBSS
- Micropipettes

Preparation of Solutions

Preparation of SPiDER-βGal DMSO stock solution

Add 14 μl of DMSO to a tube of SPiDER-βGal and dissolve it with pipetting. Store the SPiDER-βGal DMSO stock solution at -20°C.

Preparation of Bafilomycin A1 DMSO stock solution

Add 24 μl of DMSO to a tube of Bafilomycin A1 and dissolve it with pipetting. Store the Bafilomycin A1 DMSO stock solution at -20°C.

General Protocol

— Assay for living cells —

Preparation of Bafilomycin A1 working solution

Dilute the Bafilomycin A1 DMSO stock solution 1,000 times with culture medium or HBSS.

Preparation of SPiDER-βGal working solution

Mix the SPiDER-βGal DMSO stock solution and Bafilomycin A1 DMSO stock solution. Dilute the mixture 1,000 times with culture medium or HBSS.

※ For example, in order to prepare 1 ml of the SPiDER-βGal working solution, mix 1 μl of SPiDER-βGal DMSO stock solution and 1 μl of Bafilomycin A1 DMSO stock solution. Dilute the mixture with 1 ml of culture medium or HBSS.

1. Prepare cells on 35 mm dish for assay.
2. Discard the culture medium and wash the cells with 2 ml of culture medium or HBSS once.
3. Add 1 ml of Bafilomycin A1 working solution and incubate at 37°C for 1 hour in a 5%CO₂ incubator.
4. Add 1 ml of SPiDER-βGal working solution and incubate at 37°C for 30 minutes in a 5%CO₂ incubator.
5. After removing the supernatant, wash the cells with 2 ml of culture medium or HBSS twice.
6. Observe the cells under a fluorescence microscope or analyze by a flow cytometer.

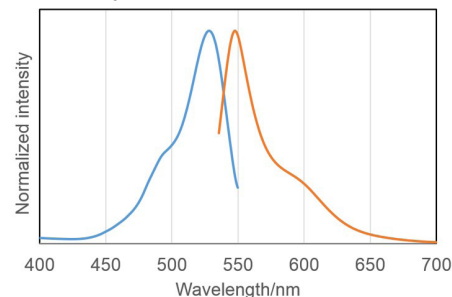


Fig. 2 Excitation and emission spectra of SPiDER-βGal after reaction with β-galactosidase

— Assay for fixed cells —

Preparation of SPiDER-βGal working solution

Dilute the SPiDER-βGal DMSO stock solution 2,000 times with Mcllvaine buffer (pH 6.0).

※ Preparation of Mcllvaine buffer (pH 6.0):

Mix 0.1 mol/l citric acid solution (3.7 ml) and 0.2 mol/l sodium phosphate solution (6.3 ml). Confirm the pH is 6.0.

If the pH is not 6.0, adjust the pH by adding either citric acid solution or sodium phosphate solution. Dilute this buffer 5 times with ultrapure water.

1. Prepare cells on 35 mm dish for assay.
2. After removing the culture medium, wash the cells with 2 ml of HBSS once. Add 2 ml of 4%paraformaldehyde (PFA) / PBS solution to the cells and incubate at room temperature for 3 minutes.
3. Remove the supernatant, and wash the cells with 2 ml of HBSS three times.
4. Add 2 ml of SPiDER-βGal working solution and incubate at 37°C for 30 minutes.

※ We recommend not to use a 5%CO₂ incubator for fixed cell experiments.

If incubation is done in a 5%CO₂ incubator, the pH of the buffer may become acidic. Acidic pH results in higher background from the endogenous β-galactosidase activity and it would be difficult to distinguish between normal cells and senescent cells.

5. After removing the supernatant, wash the cells with 2 ml of HBSS twice.
 6. Observe the cells under a fluorescence microscope or analyze by a flow cytometer.
- ※ Examples of a flow cytometry analysis are available in the product web page.

Usage Examples

Fluorescence imaging of SA-β-gal

1. WI-38 cells (5×10⁴ cells/dish, MEM, 10%fetal bovine serum, 1%penicillin-streptomycin) of passage number 0 and 12 were seeded respectively in a μ-dish 35 mm (ibidi) and cultured overnight in a 5%CO₂ incubator.
2. The cells were washed with 2 ml of HBSS once.
3. Bafilomycin A1 working solution (1 ml) was added to the culture dish, and the cells were incubated for 1 hour in a 5%CO₂ incubator.
4. SPiDER-βGal working solution (1 ml) was added to the culture dish, and the cells were incubated for 30 minutes in a 5%CO₂ incubator.
5. After the supernatant was removed, the cells were washed with 2 ml of HBSS twice.
6. HBSS (2 ml) were added and the cells were observed by confocal fluorescence microscopy (Excitation: 488 nm, Emission: 500-600 nm).

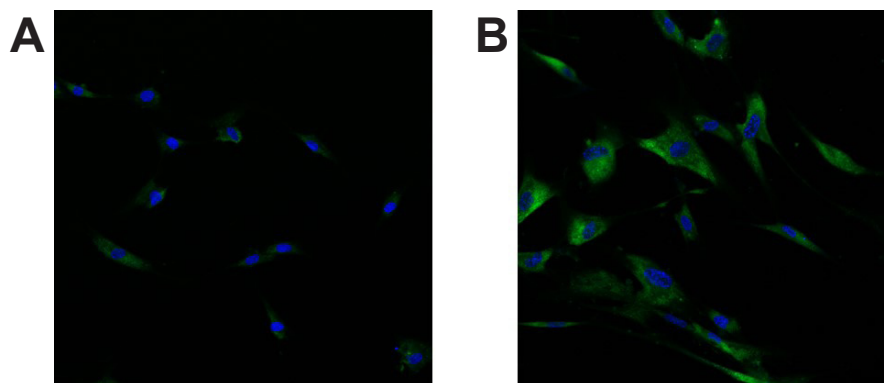


Fig. 3 Fluorescence imaging of SA-β-gal in WI-38 cells

A. Passage 0, B. Passage 12

(green: SPiDER-βGal, blue: Hoechst 33342)

Quantitative analysis of SA-β-gal positive cells by flow cytometry

1. WI-38 cells (1×10⁵ cells/dish, MEM, 10%fetal bovine serum, 1% penicillin-streptomycin) of passage number 1 and 12 were seeded respectively in a μ-dish 35 mm (ibidi) and cultured overnight in a 5%CO₂ incubator.
2. The cells were washed with 2 ml of HBSS once.
3. Bafilomycin A1 working solution (1 ml) was added to the culture dish, and the cells were incubated for 1 hour in a 5%CO₂ incubator.
4. SPiDER-βGal working solution (1 ml) was added to the culture dish, and the cells were incubated at for 30 minutes in a 5%CO₂ incubator.
5. After the supernatant was removed, the cells were washed with 2 ml of HBSS twice.
6. The cells were harvested by trypsin and resuspended in MEM (10%fetal bovine serum, 1%penicillin-streptomycin).
7. The cells were observed by a flow cytometer (Excitation: 488 nm, Emission: 515-545 nm).

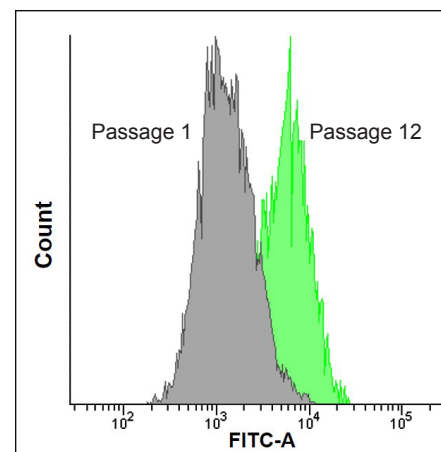


Fig. 4 Quantification of SA-β-gal positive WI-38 cells

Reference

- 1) T. Doura, M. Kamiya, F. Obata, Y. Yamaguchi, T. Y. Hiyama, T. Matsuda, A. Fukamizu, M. Noda, M. Miura and Y. Urano, *Angew. Chem. Int. Ed.*, **2016**, *55*, 9620.

For research use only. Not for use in diagnostic procedures.

If you need more information, please contact Dojindo technical service.

Dojindo Laboratories

2025-5 Tabaru, Mashiki-machi, Kamimashiki-gun, Kumamoto
861-2202, Japan Phone: +81-96-286-1515 Fax: +81-96-286-1525
E-mail: info@dojindo.co.jp Web: www.dojindo.co.jp

Dojindo Molecular Technologies, Inc.

Tel: +1-301-987-2667 Web: http://www.dojindo.com/

Dojindo EU GmbH

Tel: +49-89-3540-4805 Web: http://www.dojindo.eu.com/

Dojindo China Co., Ltd

Tel: +86-21-6427-2302 Web: http://www.dojindo.cn/